



The easy way to secure vessels on shaker platforms without the use of clamps




Thermo Scientific Adhesive Flask Mat

Secure your vessels with confidence

- **Easy-to-view samples**
The white background makes it easy to determine when to harvest your samples
- **Simple to clean**
Can be restored to like new condition by wiping with alcohol
- **Long lasting**
Months of use without losing adhesive properties

Our Adhesive Flask Mat is an alternative to flask clamps for securing odd shaped objects such as bottles, volumetric flasks and test tube racks to a shaker platform. Simply place the 9" x 9" Adhesive Flask Mat on an anodized aluminum or stainless steel shaker platform and secure the vessel to it. It will hold Erlenmeyer flasks as large as 4 liters that contain 2 liters of media at speeds up to 250 rpm. The mat is not suitable for use in water baths.

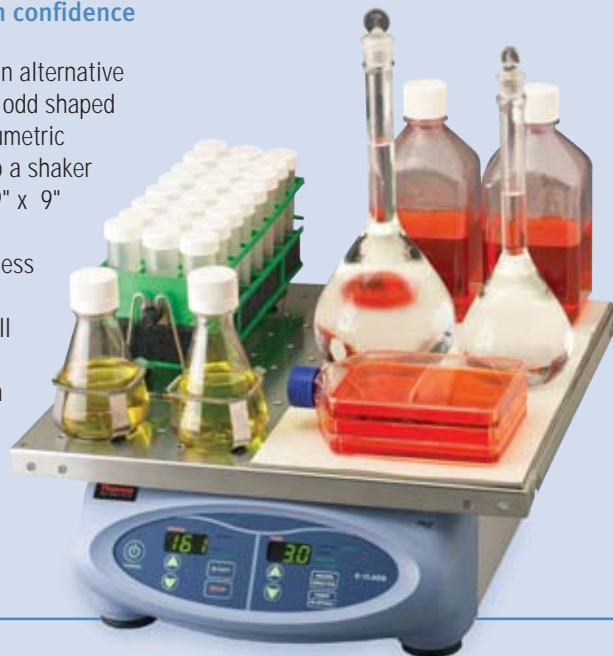

Adhesive Flask Mat

Cat. No.	Description	Quantity
300349	9" x 9" Adhesive Flask Mat	each



North America: USA/Canada +1 866 984 3766

Europe: Austria +43 1 801 40 0, Belgium +32 2 482 30 30, Finland +358 9 329 100, France +33 2 2803 2000, Germany national toll free 08001-536 376, Germany international +49 6184 90 6940, Italy +39 02 95059 1, Netherlands +31 76 571 4440, Russia/CIS +7 (812) 703 42 15, Spain/Portugal +34 93 223 09 18, Switzerland +41 44 454 12 12, UK/Ireland +44 870 609 9203

Asia: China +86 21 6865 4588 or +86 10 5850 3588, India +91 22 5542 9494, Japan +81 45 453 9220, Other Asian countries +852 2885 4613

Countries not listed: +49 6184 90 6940 or +33 2 2803 2000

www.thermo.com/shakers

DS-LESHKFLMAT-0408

© 2008 Thermo Fisher Scientific Inc. All rights reserved. All trademarks are the property of Thermo Fisher Scientific Inc. and its subsidiaries. Specifications, terms and pricing are subject to change. Not all products are available in all countries. Please consult your local sales representative for details.

Thermo
SCIENTIFIC